

ITGB7 Monoclonal Antibody

Description

Product type	Antibody
Code	BT-MCA3092
Host	Mouse
Isotype	Mouse IgG1
Size	100µL, 50µL
Immunogen	Purified recombinant fragment of human ITGB7 (AA: extra 20-174) expressed in E. Coli.
Mol wt	86.9kDa
Species reactivity	Human
Clonality	Monoclonal
Recommended application	IHC,FCM
Concentration	N/A
Full name	N/A
Synonyms	N/A

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

Background

This gene encodes a protein that is a member of the integrin superfamily. Members of this family are adhesion receptors that function in signaling from the extracellular matrix to the cell. Integrins are heterodimeric integral membrane proteins composed of an alpha chain and a beta chain. The encoded protein forms dimers with an alpha4 chain or an alphaE chain and plays a role in leukocyte adhesion. Dimerization with alpha4 forms a homing receptor for migration of lymphocytes to the intestinal mucosa and Peyer's patches. Dimerization with alphaE permits binding to the ligand epithelial cadherin, a calcium-dependent adhesion molecule. Alternate splicing results in multiple transcript variants. Additional alternatively spliced transcript variants of this gene have been described, but their full-length nature is not known.

Recommended Dilution

WB: 1:500 - 1:2000


IHC-p: 1:200 - 1:1000

FCM: 1:200 - 1:400


ELISA: 1:10000

Not yet tested in other applications.


Images

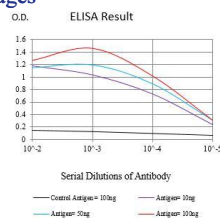


Black line: Control Antigen (100 ng);Purple line: Antigen (10ng); Blue line: Antigen (50 ng); Red line:Antigen (100 ng)

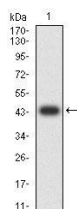


Western blot analysis using ITGB7 mAb against human ITGB7 (AA: extra 20-174) recombinant protein. (Expected MW is 43.4 kDa)

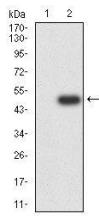


Western blot analysis using ITGB7 mAb against HEK293-6e (1) and ITGB7 (AA: extra 20-174)-hlgGfc transfected HEK293-6e (2) cell lysate.

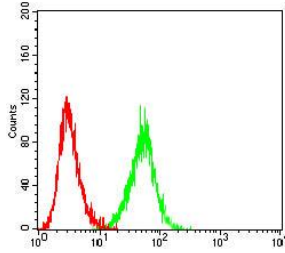

Flow cytometric analysis of Jurkat cells using ITGB7 mouse mAb (green) and negative control (red).

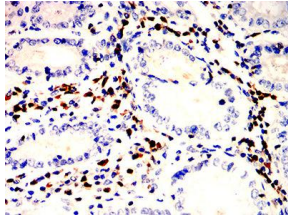

Immunohistochemical analysis of paraffin-embedded gastric cancer tissues using ITGB7 mouse mAb with DAB staining.

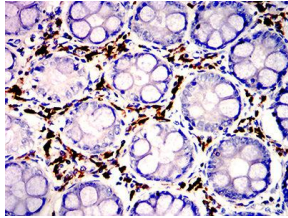

Immunohistochemical analysis of paraffin-embedded colon tissues using ITGB7 mouse mAb with DAB staining.

Storage

Store at 4°C short term. Aliquot and store at -20°C long term.

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: save@bt-laboratory.com | www.bt-laboratory.com