

## CHRNA7 Monoclonal Antibody

### Description

<b>Product type</b>	Antibody
<b>Code</b>	BT-MCA4265
<b>Host</b>	Mouse
<b>Isotype</b>	Mouse IgG1
<b>Size</b>	100µL, 50µL
<b>Immunogen</b>	Purified recombinant fragment of human CHRNA7 (AA: extra 52-259) expressed in E. Coli.
<b>Mol wt</b>	56.4kDa
<b>Species reactivity</b>	Human
<b>Clonality</b>	Monoclonal
<b>Recommended application</b>	IHC,FCM
<b>Concentration</b>	N/A
<b>Full name</b>	N/A
<b>Synonyms</b>	NACHRA7;CHRNA7-2

This product is for research use only, not for use in human, therapeutic or diagnostic procedure.

### Background

N/A

### Recommended Dilution

WB: 1:500 - 1:2000


IHC-p: 1:200 - 1:1000

FCM: 1:200 - 1:400


ELISA: 1:10000

Not yet tested in other applications.


### Images

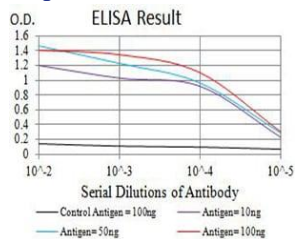


Black line: Control Antigen (100 ng);Purple line: Antigen (10ng); Blue line: Antigen (50 ng); Red line:Antigen (100 ng)

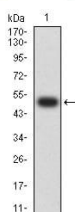


Western blot analysis using CHRNA7 mAb against human CHRNA7 (AA: extra 52-259) recombinant protein. (Expected MW is 50.4 kDa)

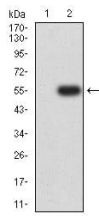


Western blot analysis using CHRNA7 mAb against HEK293 (1) and CHRNA7 (AA: extra 52-259)-hlgGfc transfected HEK293 (2) cell lysate.

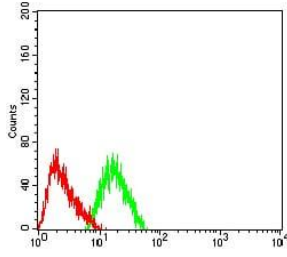

Flow cytometric analysis of SH-SY5Y cells using CHRNA7 mouse mAb (green) and negative control (red).

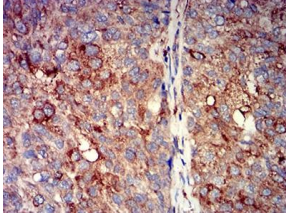

Immunohistochemical analysis of paraffin-embedded bladder cancer tissues using CHRNA7 mouse mAb with DAB staining.

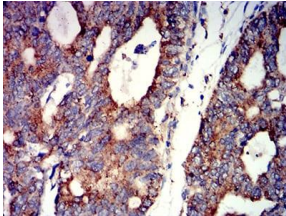

Immunohistochemical analysis of paraffin-embedded rectum cancer tissues using CHRNA7 mouse mAb with DAB staining.

### Storage

Store at 4°C short term. Aliquot and store at -20°C long term.

501 Changsheng S Rd, Nanhu Dist, Jiaxing, Zhejiang, China

Tel: 86 21 31007137 | E-mail: [save@bt-laboratory.com](mailto:save@bt-laboratory.com) | [www.bt-laboratory.com](http://www.bt-laboratory.com)